

# Inno-Hook™

## Polymer Ligation System

- Compatible with the corresponding sizes of Inno-Hook polymer ligating clips

**Vessel/Tissue Bundle Size Range**

Polymer Ligating Clips, ML			Medium Large for 3 mm-10 mm		5mm, ML
Polymer Ligating Clips, L			Large for 5 mm-13 mm		10mm, L
Polymer Ligating Clips, XL			Extra Large for 7 mm-16 mm		10mm, XL


*Reusable Clip Applicator (Detachable), Single-Fire*




- Wider clip aperture enables ligation of larger structure than insertion port



- Unique detachable forehead reduces the maintenance cost due to damaged jaws

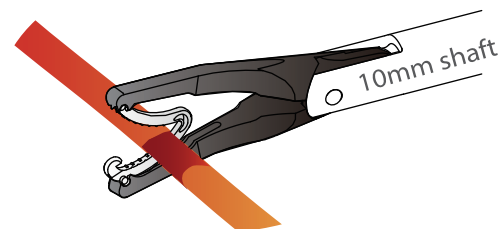

up to 16mm Vessels

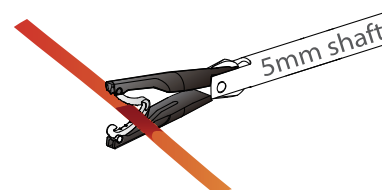

up to 10mm Vessels

- Color-coded knobs matches with the cartridge tints

- Stainless steel with specific metal process for durability

- Flushing port facilitates cleaning process to prevent potential contamination



- Ergonomic design allows one-handed rotation and operation

*Reusable Clip Applicator (Detachable), Single-Fire, 5mm, ML (#RD5P1-B0)*

Security on Structure

Inno-Hook™ polymer ligating clips are designed with the locking mechanism to fully close the structure and securely withstand clip slipping and burst pressure force.

Demonstrate better retention performance with unique inner texture against transverse and axial force by the empirical comparison with Hem-o-lok.

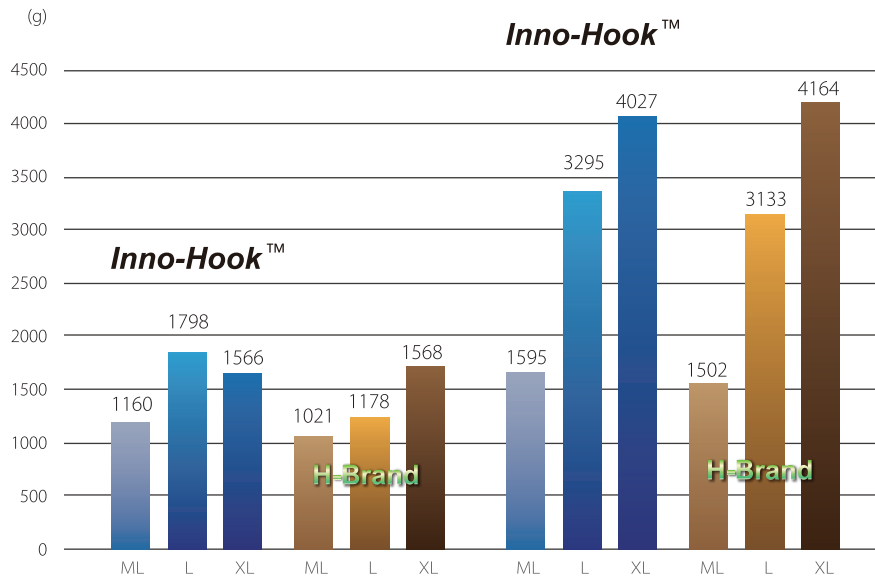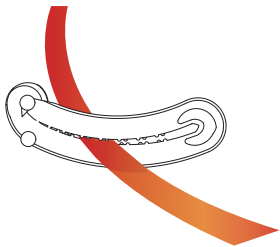

Safety in Essence

Inno-Hook™ non-absorbable polymer clips pass biocompatibility tests in accordance with ISO 10993-1 to ensure the safety for human implantation. Immune from the thermal damages risk caused by the energy-based ligation products.

Inno-Hook™ Reusable Clip Applier

Product No.	Description	Color	Shaft Spec. (D*U)	Packaging
RA5P1-B0	Reusable Clip Applier (Standard), Single-Fire, 5mm, Polymer, ML	Green	ø5*320mm	1 pcs / box, 5 boxes / carton
RAXP1-L0	Reusable Clip Applier (Standard), Single-Fire, 10mm, Polymer, L	Purple	ø10*320mm	1 pcs / box, 5 boxes / carton
RAXP1-X0	Reusable Clip Applier (Standard), Single-Fire, 10mm, Polymer, XL	Brown	ø10*320mm	1 pcs / box, 5 boxes / carton
RD5P1-B0	Reusable Clip Applier (Detachable), Single-Fire, 5mm, Polymer, ML	Green	ø5*320mm	1 pcs / box, 5 boxes / carton
RDXP1-L0	Reusable Clip Applier (Detachable), Single-Fire, 10mm, Polymer, L	Purple	ø10*320mm	1 pcs / box, 5 boxes / carton
RDXP1-X0	Reusable Clip Applier (Detachable), Single-Fire, 10mm, Polymer, XL	Brown	ø10*320mm	1 pcs / box, 5 boxes / carton
Product No.	Description	Color	SHAFT Spec. (Diameter)	Packaging
AD5P1-B0	Detachable Jaw, 5mm, Polymer Clip, ML	Green	ø5mm	4 pcs / box
ADXP1-L0	Detachable Jaw, 10mm, Polymer Clip, L	Purple	ø10mm	4 pcs / box
ADXP1-X0	Detachable Jaw, 10mm, Polymer Clip, XL	Brown	ø10mm	4 pcs / box

Inno-Hook™ Polymer Clips

Product No.	Description	Color	Vessel/Tissue Size Range	Packaging
AA5P1-BC	Non-absorbable Polymer Ligating Clips, ML, 3 clips	Green	3mm-10mm	14 cartridges / box, 12 boxes / carton
AA5P1-BF	Non-absorbable Polymer Ligating Clips, ML, 6 clips	Green	3mm-10mm	14 cartridges / box, 12 boxes / carton
AAXP1-LC	Non-absorbable Polymer Ligating Clips, L, 3 clips	Purple	5mm-13mm	14 cartridges / box, 12 boxes / carton
AAXP1-LF	Non-absorbable Polymer Ligating Clips, L, 6 clips	Purple	5mm-13mm	14 cartridges / box, 12 boxes / carton
AAXP1-XC	Non-absorbable Polymer Ligating Clips, XL, 3 clips	Brown	7mm-16mm	14 cartridges / box, 12 boxes / carton
AAXP1-XF	Non-absorbable Polymer Ligating Clips, XL, 6 clips	Brown	7mm-16mm	14 cartridges / box, 12 boxes / carton

